

# SPECIMEN COLLECTION CHART

**All specimens must be labeled with a first and last name and a second identifier, e.g., date of birth.**




## Pap/HPV Testing: ThinPrep

- HPV Only
- HPV Only Reflex to Genotype
- Pap 1: Pap Test
- Pap 2: Pap Test Reflex to HPV if ASCUS
- Pap 3: Pap Test and HPV
- Pap 4: Pap Test w/ HPV, Reflex to Genotypes 16, 18/4
- Pap 5: Pap Test Reflex to HPV if ASCUS or AGUS
- Pap 6: Pap Test Reflex to HPV if ASCUS or AGUS, Reflex to Genotype 16, 18/45
- Pap 7: Pap Test Reflex to HPV if ASCUS, ASC-H, or AGUS
- Pap 8: Pap Test Reflex to HPV if ASCUS, ASC-H, or AGUS, Reflex to Genotype 16, 18/45
- Pap 9: Pap Test Reflex to HPV if ASCUS, Reflex to Genotype 16, 18/45
- Pap 10: Pap Test, ThinPrep, Image Guided, Reflex to HPV if ASCUS, ASC-H, AGUS, or LSIL
- Pap 11: Pap Test, ThinPrep, Image Guided, Reflex to HPV if ASCUS, ASC-H, AGUS, or LSIL, Reflex Genotype 16, 18/45

**Room Temp**

## Myswab

- Bacterial Vaginosis - Swab
- Candida Vaginitis - Swab
- Chlamydia trachomatis - Swab
- Herpes Simplex Virus 1 and 2 (HSV) - Swab
- MySwab STD Profile - Swab
- MySwab Vaginosis Profile + STD Profile - Swab
- MySwab Vaginosis Profile + Trichomonas - Swab
- MySwab Vaginosis Profile - Swab
- Neisseria gonorrhoeae - Swab
- Trichomonas - Swab

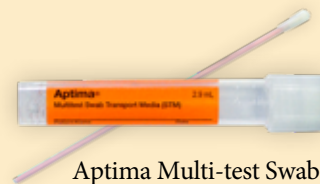

Aptima Multi-test Swab

\*5 tests performed per swab \*

**Room Temp**

## Vaginal/Genital Cultures:

### Genital/Vaginal Cultures, Group B Streptococcus

#### Specimen Options:

- |   |  |
|---|--|
| White Copan Eswab<br>(Preferred)<br><span style="color: red;">Break swab inside of container and replace cap</span> | Green Copan Eswab<br>"Flexible minitip" (Alternative)<br><span style="color: red;">Break swab inside of container and replace cap</span> |
|---|--|



White Copan Eswab

Green Copan Eswab

**Room Temp**

## Vaginitis Testing:



BD Affirm Ambient Temperature Transport Swabs (ATTS)

### Testing includes:

- Gardnerella
- Candida (7 species)
- Trichomonas

**Room Temp**

## Vaginal Profile Swab, no pH



White Copan Eswab (Preferred)

Green Copan Eswab  
"Flexible minitip" (Alternative)

### Testing includes:

- Trichomonas Antigen
- Wet Prep for Yeast
- Whiff Test
- Gram Stain for Bacterial Vaginosis
- Interpretation

**Room Temp**

## Vaginal Profile Saline, with pH, Outreach



Sterile container with 1.5 mLs of saline


### Testing includes:

- Trichomonas Antigen
- Wet Prep for Yeast
- pH
- Whiff Test
- Gram Stain for Bacterial Vaginosis
- Interpretation

**Room Temp**

## STI Testing by Molecular: Chlamydia, Neisseria Gonorrhoeae, and Trichomonas: by Molecular

### Specimen Options:

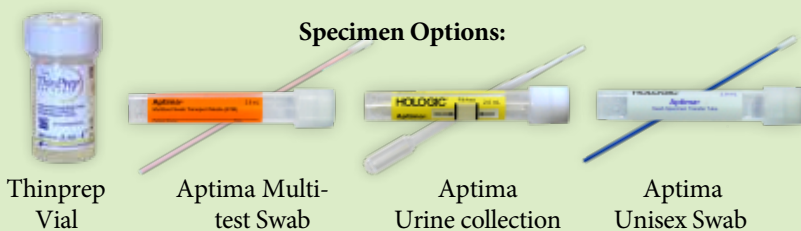

ThinPrep Vial

Aptima Multi-test Swab

Aptima Urine collection

Aptima Unisex Swab

**Room Temp**

## STI- Herpes Simplex Virus 1 and 2 (HSV) Swab - TMA DNA, PCR



UTM Mycoplasma genitalium

**Refrigerated**



Aptima Multi-test Swab

**Room Temp**

